



**Test system for function tests of clinical thermometer
(DIN EN ISO 80601-2-56), cardiac output-, invasive (IEC 60601-2-34)
and non-invasive blood pressure devices (IEC 80601-2-30)**



- ☑ temperature simulation from 0 – 70 °C
- ☑ simulation of multiple NIBP signals
- ☑ suitable for calibration of medical thermometer
- ☑ user specific language setting
- ☑ accumulator operated device in a carrying case
- ☑ dynamic and static IBP signal simulation
- ☑ 4 different cardiac outputs selectable
- ☑ Versions available: NE, NT, E, N

Technical Data

Mains voltage:	83 – 264 V ac, 50 / 60 Hz or internal accumulator operation	Static signals:	0 - 400 mmHg in 1 mmHg steps
Nominal power:	max. 25 VA	Dynamic signals:	± 1 mmHg
Protection class:	internal power supply	ARP 120 / 80 mmHg	aorta radial pressure
Environmental temperature:	+ 5 - + 40 °C	AFP 130 / 70 mmHg	aorta femoral pressure
Storage temperature:	- 10 - + 50 °C	ZVP 10 / 0 mmHg	central vein pressure
Timer:	10 – 60 min in 10 min steps, ± 1 sec	PAP 25 / 10 mmHg	pulmonary aorta pressure
Interface:	1 x RS-232 for PC connection	PAWP 10 / 0 mmHg	pulmonary-ast-wedge-position
Digital display:	4 x 16 char display	LVP 10 / 0 mmHg	left ventricle pressure
Keyboard:	8 key foil keyboard	LAP 5 / 0 mmHg	left atrium pressure
Accessories:	1 x RS-232 interface cable, charger or power supply	RVP 25 / 0 mmHg	right ventricle pressure
Mechanical data:	light weight metal case IP20	RAP 5 / 0 mmHg	right atrium pressure
Dimensions:	138 x 230 x 37 mm (W x H x D)	Arrhythmia:	aorta damped 110 / 80 mmHg
Weight:	approx. 6,5 kg (incl. suitcase and accessories)	Puls rates:	aorta disturbed 170 / 20 mmHg
			30 - 150 bpm in 10 bpm steps
			± 2 bpm
Temperature:	0 - 70 °C in 1 °C steps, ± 0,03 °C	Non-invasive blood pressure:	0,1 mmHg steps
	20 - 50 °C in 0,1 °C steps, ± 0,03 °C	Manometer function:	0 - 300 mmHg, ± 0.8 mmHg
	YSI 400 / 700 compatible	Leakage rate measurement:	0 - 300 mmHg, ± 0.8 mmHg
Temperature [°C]:	20,40 / 21,70 / 23,00 / 24,30 / 25,60 / 26,90 28,20 / 29,50 / 30,80 / 32,10 / 33,40 / 34,70 36,00 / 37,30 / 38,60 / 39,90 / 41,20 / 42,50 43,80 / 45,10 / 46,40 / 47,70 / 49,00 ± 0,03 °C	Simulation mode:	1 - 15 min in 1 min steps, ± 1 sec
Cardiac output:	3, 5, 7, 9 l/min, ± 0,3 l/min	Puls rate:	oscillometric
Blood temperature:	37 °C, ± 0,1 °C	Reference volume:	80 bpm, ± 2 bpm
Injectate temperature:	1, 4, 12, 20 °C, ± 0,2 °C	Systole / Diastole:	100 / 500 ml
Injectate volume:	10 ml	Time Measurement	120/80, 180/120, 90/50 mmHg ± 3 mmHg
		Accessories:	1 – 1000 sec, ± 0,1 sec, ± 1 %
			hand calibration pump
			adapter block with connection set
			reference pressure volume
			100 / 400 ml
Invasive blood pressure:	5 or 40 µV/ V/ mmHg [V(monitor)]	Selectable languages:	german, english, french, polish
Output impedance:	178, 345 Ohm, ± 2 Ohm		spanish, italian, portuguese, turkish

Temperature – simulator

The simulation of temperature levels is based on precision resistances with fixed output values, that can be adjusted in 1°C steps (0-70°C) and 0,1°C steps (20-50°C). Accurate resistance values, like for the thermistorprobe series YSI 400/ 700, are implemented for simulation.

Cardiac output – simulator

By changing the temperature / resistance ratio send to the monitor, the cardiac output can be simulated. You can choose between four different output values and four injectate temperatures. The blood temperature is fixed at 37°C and the injectate volume to 10ml.

Invasive – blood pressure – simulator


Static and dynamic signals can be simulated by changing the values send from the IBP modules to the transducer. You can choose between various signal styles and arrhythmia. In addition to that, the pulse rate can be adjusted in 10bpm steps from 30 to 150 bpm.

Non-invasive – blood pressure – simulator

The blood pressure simulator comes with a precision hand pump, an adapter block, hoses and two reference pressure volumes. Static pressures can be generated via the hand pump and are simultaneously displayed in the simulator. Three different systole / diastole combinations can be simulated via an oscillometric method. The pulse rate is fixed at 80bpm. There is a leakage rate measurement and a manometer function included.

(The specified measuring accuracy refers to the measuring element. Technical modifications and errors reserved. 11/2022)

Technical Data

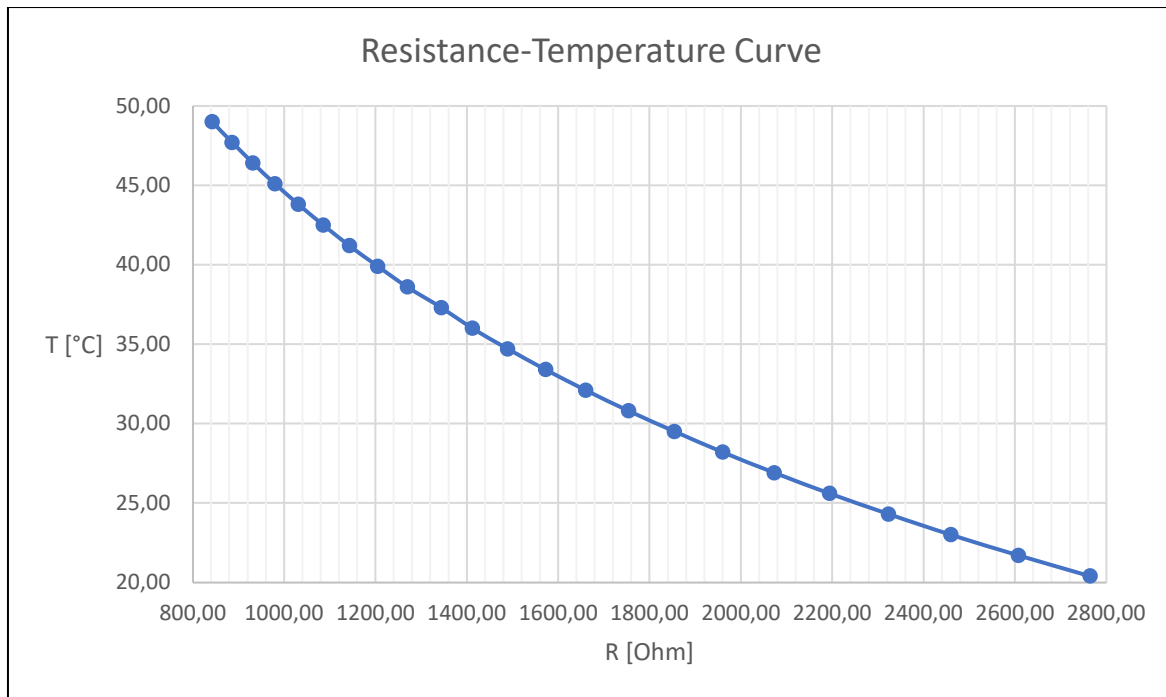

R [Ohm]	T [°C]
2764,40	20,40
2607,40	21,70
2460,00	23,00
2323,40	24,30
2194,40	25,60
2073,20	26,90
1960,40	28,20
1854,50	29,50


R [Ohm]	T [°C]
1754,20	30,80
1660,20	32,10
1572,60	33,40
1489,60	34,70
1412,00	36,00
1344,20	37,30
1269,80	38,60
1204,90	39,90

R [Ohm]	T [°C]
1143,00	41,20
1085,50	42,50
1031,20	43,80
979,59	45,10
931,80	46,40
886,12	47,70
842,80	49,00


R [Ohm] / T [°C]
 6989 / 1,0
 3886 / 13,0
 2252 / 25,0
 1355 / 37,0
 842,8 / 49,0

(The specified measuring accuracy refers to the measuring element. Technical modifications and errors reserved. 11/2022)


Technical Data

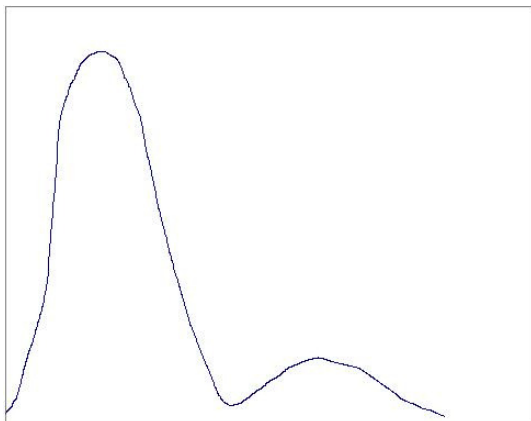

aortafemoral pressure
130 / 70 mmHg (AFP)

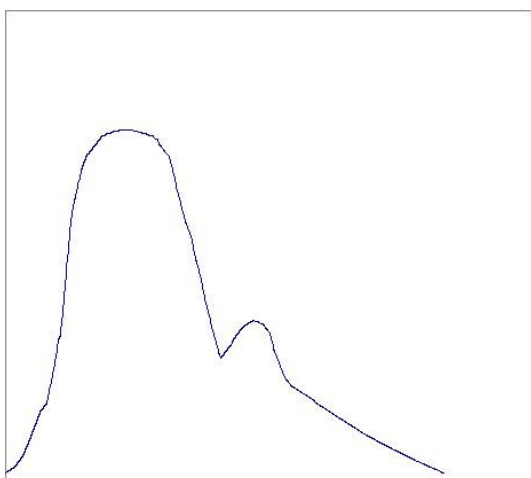


aortaradial pressure
120 / 70 mmHg (ARP)

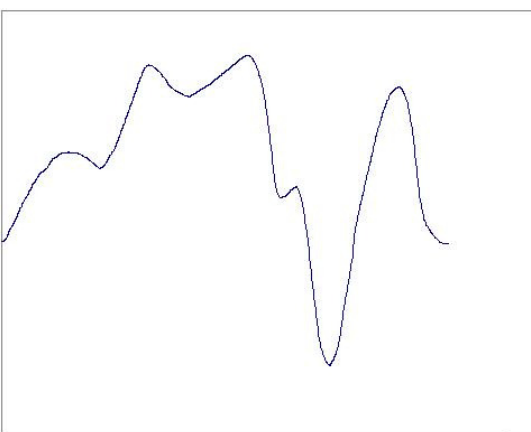


central venous pressure
10 / 0 mmHg (ZVP)

(The specified measuring accuracy refers to the measuring element. Technical modifications and errors reserved. 11/2022)


Technical Data

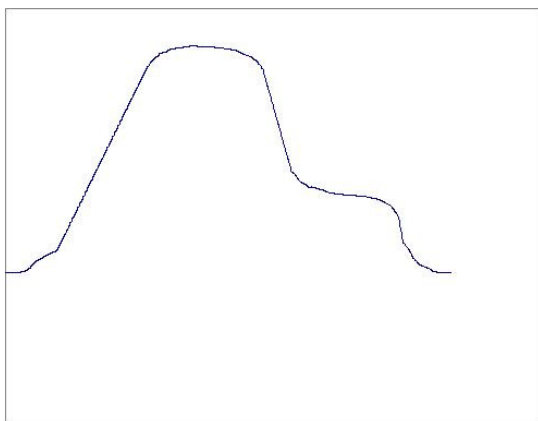

pulmonary pressure
25 / 10 mmHg (PAP)

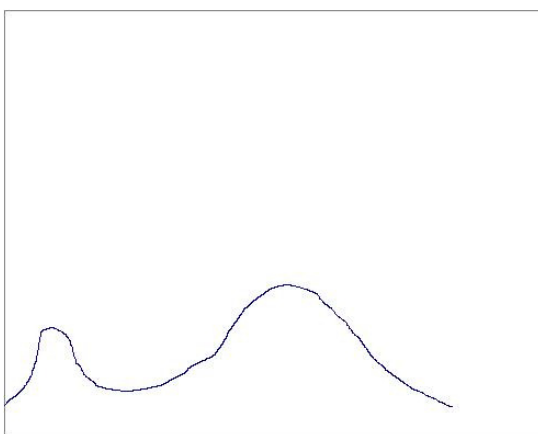


pulmonalis branch, wedge-position
10 / 0 mmHg (PAWP)

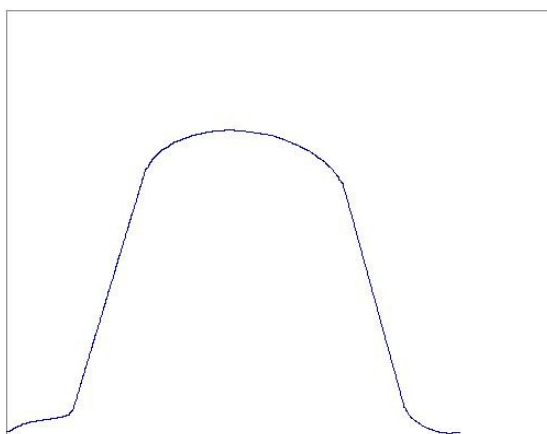


left ventricle pressure
10 / 0 mmHg (LVP)

(The specified measuring accuracy refers to the measuring element. Technical modifications and errors reserved. 11/2022)


Technical Data

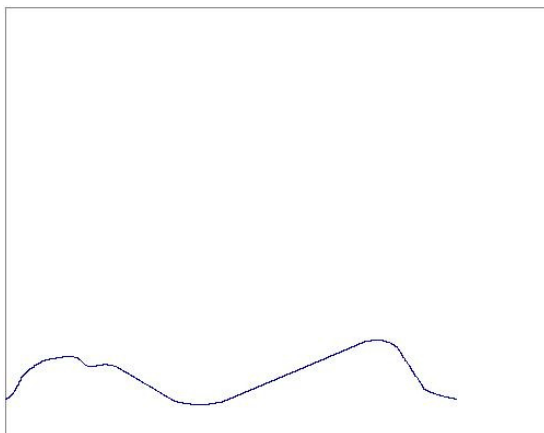

left intraarterial pressure
5 / 0 mmHg (LAP)

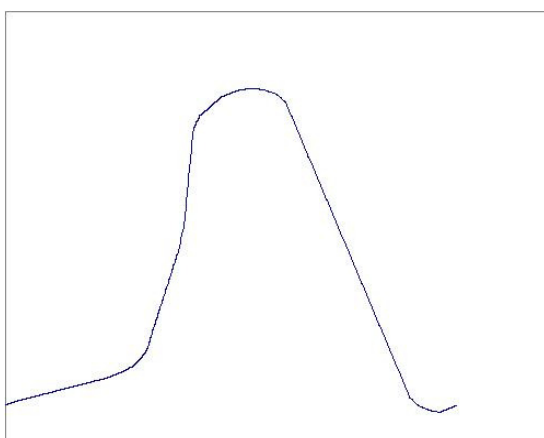

right intraventricular pressure
25 / 0 mmHg (RVP)

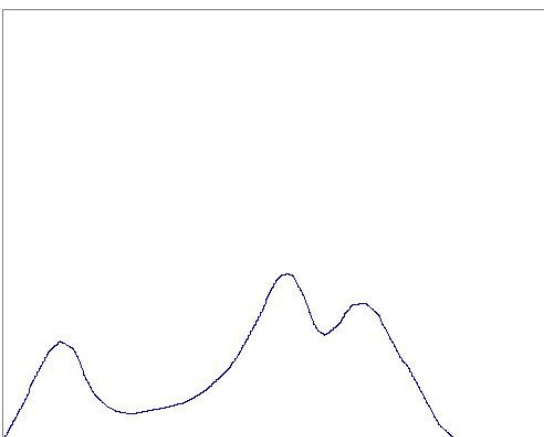

right intraarterial pressure
5 / 0 mmHg (RAP)

(The specified measuring accuracy refers to the measuring element. Technical modifications and errors reserved. 11/2022)

Example of use

Versions:

	Version			
	NE	NT	E	N
T NIBP IBP CO	☑	☑	☑	
	☑	☑		☑
	☑		☑	
	☑		☑	

T: Temperature
NIBP: Non-invasive blood pressure
IBP: Invasive blood pressure
CO: Cardiac Output

(The specified measuring accuracy refers to the measuring element. Technical modifications and errors reserved. 11/2022)